



We put space to work.

CaseVue™ Preconfigured Carts

One look and the advantage is clear.

Built for easy transport, the CaseVue™ case cart is stronger, lighter, and the most ergonomic choice available. Carts can be used wherever contents must be protected from an “unclean” environment, including hospital based or offsite surgery centers.

Advantages provided by the patent pending CaseVue case carts include:

Proficient: Maintenance-free casters reduce the total cost of ownership. Optional clear doors and status flags provide a quick view of cart's status at any time. Cool touch handles eliminate the need for protective cloth to remove cart from cart wash. Biohazard wheel indicates used carts and paper holder conveniently holds pick lists and routing slips.

Clean Design: Cart is designed to be easily cleaned and not retain water. Stainless steel corner handles are encased in a cool touch polymer featuring Microban® antimicrobial product protection.

Less Effort: Premier, quiet, maintenance-free, 6" diameter casters reduce rolling resistance; wide grip handles provide better leverage for steering and swivel lock casters keep the cart on track. Optional upgrades, including 5th wheel steering and swing up handles, provide maximum control in transit.

Quiet: Contoured sides, reinforced top and floor, softer wheels and other noise dampeners, provide a noticeably quieter ride.

Durable: Type 304 stainless steel construction, including stainless steel casters are standard, not an optional upgrade.

Efficient: Metro CaseVue delivers the most storage per square foot to ensure we are making the most of your space. Roller shelves extend further to provide better access.

Configurable: CaseVue provides many standard configurations to select the right cart for your needs, including clear doors, various sizes, solid or wire roller shelves, 5th wheel steering, colored status flags and swing up handle.



CV24HC-NBFH-S2



CV36LS-5BFH-W1



CaseVue™ Preconfigured Carts



LO2-251
Printed in U.S.A. 1/22

Information and specifications are subject to change without notice. Please confirm at time of order.

Copyright © 2022 InterMetro Industries Corp.

33.12



Specifications

Cart Body

- Type 304 Stainless Steel, No. 3 finish type
- Contoured side panels to rigidize body and dampen vibration noise
- Four type 304 stainless steel, wide grip handles encased in injection molded flexible PVC with Microban antimicrobial product protection
- Sloped flanges to promote water drainage
- Raised front on floor to retain spills in transit
- Reinforced top and floor to stiffen surface and reduce vibration noise
- Stainless steel door stay with impact-resistant polymer catch
- Polymer and stainless steel roll pin paper holder

Doors

- Type 304 stainless steel single panel
- Optional: polycarbonate clear window sealed with rubber gasket
- Sloped flanges to promote drainage
- Replaceable paddle latch with security option
- Rubber bumpers to reduce rattle noise

Base Assembly

- Replaceable non marking, type 304 stainless corner bumpers encased in flexible PVC, with front and rear continuous bumper
- Replaceable type 304 stainless steel caster channels

Shelves

- Roller pullout stainless steel wire or solid
- Built in retainer that holds shelf closed during cart wash
- Each shelf includes two shelf brackets with nylon rollers
- Shelf bracket include a shelf stop to prevent over extension
- Adjustable on 2.25" (57mm) increments
- Rubber bumpers to minimize rattle noise

Swing Up Extended Handle

- Type 304 stainless steel tubular construction
- Swings up to desired height of user
- Closes within footprint of the cart
- Rubber bumpers reduce rattle noise

Available options

(Minimums may apply)

- Extended tops
- Tote boxes
- Custom sizes
- Tow assemblies

Casters

Cart washable, Quiet, Easy Rolling, Maintenance-Free

- 6" (152mm) diameter, type 304 stainless steel
- Two total-lock casters on front of cart
- Two swivel-lock on rear for carts without 5th wheel steering
- Two swivel on rear for carts equipped with 5th wheel steering
- All stainless steel hardware and top plate
- Swivel bearing: Double race of stainless steel ball bearings, grease sealed, dust cap.
- Wheel bearing: Sealed, precision, ball bearing.
- Gray thermoplastic rubber (TPR) on gray polypropylene hub with thread guard; 1¼" (32mm) rounded tread face provides a softer, quieter ride. TPR resists flat spotting and repels floor debris.

5th Wheel Steering

- Available on select carts
- 5" (127mm) diameter, type 304 stainless steel
- Type 304 stainless steel mechanism with vinyl dipped foot activated lever

Biohazard Wheel

- Type 304 stainless steel housing and wheel
- Wheel rotates and locks into position indicating Biohazard in red, or Clean in green

Color Status Flags

- 4 Injection molded ABS flags in red, green, orange and blue
- Stainless steel housing
- Flags positively held in up or down position
- Default color code (user definable)
 - Green: Cart fully picked and ready
 - Orange: Cart picked, OR supply needed
 - Red: Biohazard, soiled contents inside
 - Blue: Clean cart return (i.e. case cancelled)

Metro Closed Case Carts are cart washable. The cart wash chamber must be non-pressurized. Water temperatures in the range of 140°-180°F (60°-82°C) are allowed with an exposure period not to exceed 30 minutes. Hot air drying is allowed with the temperatures not to exceed 200°F (93° C) and an exposure period not to exceed 15 minutes. Contact your Metro representative if cart wash parameters fall outside of above specifications.




CaseVue™ Preconfigured Carts

Case Cart with Factory Installed Accessories

CaseVue standard features include, all stainless steel construction, cool touch antimicrobial wide grip handles, casters, security door latch, picklist/paper holder, door catch and optional 5th wheel steering (available on 36 & 48 models).

Model Number Identifier:

CV36LC-NBFH-S1

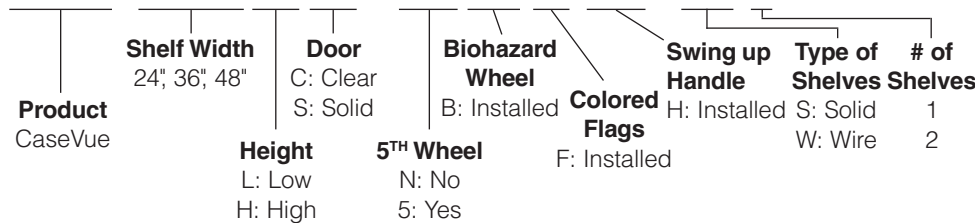

Models including **B**, **F**, or **H** have corresponding accessories factory assembled and ready to go out of the box.
B: Biohazard indicator wheel
F: Color flags
H: Swing up handle.

Preconfigured CaseVue carts include wire or solid shelves, 1 roller shelf with low carts or 2 roller shelves with tall carts. Shelves are packaged separately to be installed onsite.

Clear Door Carts

Preconfigured Carts with Shelves, Accessories and 5th Wheel

Wire Shelf	Solid Shelf	Description	D x W x H	D x W x H
CV36LC-5BFH-W1	CV36LC-5BFH-S1	Case Cart, Low, Clear Doors, 5th Wheel, Accessories, 1 Shelf	28 x 45½ x 40%	711 x 1156 x 1032
CV48LC-5BFH-W1	CV48LC-5BFH-S1	Case Cart, Low, Clear Doors, 5th Wheel, Accessories, 1 Shelf	28 x 57½ x 40%	711 x 1461 x 1032
CV36HC-5BFH-W2	CV36HC-5BFH-S2	Case Cart, High, Clear Doors, 5th Wheel, Accessories, 2 Shelves	28 x 45½ x 58¼	711 x 1156 x 1480

Preconfigured Carts with Shelves and Accessories

Wire Shelf	Solid Shelf	Description	D x W x H	D x W x H
CV24LC-NBF-W1	CV24LC-NBF-S1	Case Cart, Low, Clear Door, Accessories, 1 Shelf	28 x 33½ x 40%	711 x 851 x 1032
CV36LC-NBFH-W1	CV36LC-NBFH-S1	Case Cart, Low, Clear Doors, Accessories, 1 Shelf	28 x 45½ x 40%	711 x 1156 x 1032
CV48LC-NBFH-W1	CV48LC-NBFH-S1	Case Cart, Low, Clear Doors, Accessories, 1 Shelf	28 x 57½ x 40%	711 x 1461 x 1032
CV24HC-NBFH-W2	CV24HC-NBFH-S2	Case Cart, High, Clear Door, Accessories, 2 Shelves	28 x 33½ x 58¼	711 x 851 x 1480
CV36HC-NBFH-W2	CV36HC-NBFH-S2	Case Cart, High, Clear Doors, Accessories, 2 Shelves	28 x 45½ x 58¼	711 x 1156 x 1480

Solid Door Carts

Preconfigured Carts with Shelves, Accessories and 5th Wheel

Wire Shelf	Solid Shelf	Description	D x W x H	D x W x H
CV36LS-5BFH-W1	CV36LS-5BFH-S1	Case Cart, Low, Solid Doors, 5th Wheel, Accessories, 1 Shelf	28 x 45½ x 40%	711 x 1156 x 1032
CV48LS-5BFH-W1	CV48LS-5BFH-S1	Case Cart, Low, Solid Doors, 5th Wheel, Accessories, 1 Shelf	28 x 57½ x 40%	711 x 1461 x 1032
CV36HS-5BFH-W2	CV36HS-5BFH-S2	Case Cart, High, Solid Doors, 5th Wheel, Accessories, 2 Shelves	28 x 45½ x 58¼	711 x 1156 x 1480

Preconfigured Carts with Shelves and Accessories

Wire Shelf	Solid Shelf	Description	D x W x H	D x W x H
CV24LS-NBF-W1	CV24LS-NBF-S1	Case Cart, Low, Solid Door, Accessories, 1 Shelf	28 x 33½ x 40%	711 x 851 x 1032
CV36LS-NBFH-W1	CV36LS-NBFH-S1	Case Cart, Low, Solid Doors, Accessories, 1 Shelf	28 x 45½ x 40%	711 x 1156 x 1032
CV48LS-NBFH-W1	CV48LS-NBFH-S1	Case Cart, Low, Solid Doors, Accessories, 1 Shelf	28 x 57½ x 40%	711 x 1461 x 1032
CV24HS-NBFH-W2	CV24HS-NBFH-S2	Case Cart, High, Solid Door, Accessories, 2 Shelves	28 x 33½ x 58¼	711 x 851 x 1480
CV36HS-NBFH-W2	CV36HS-NBFH-S2	Case Cart, High, Solid Doors, Accessories, 2 Shelves	28 x 45½ x 58¼	711 x 1156 x 1480



Handle swings up to ideal height.

CaseVue™ Preconfigured Carts

We put space to work.

Model 24 (33½" L)



CV24LC-NBF-W1



CV24LS-NBF-S1



CV24HC-NBFH-S2



CV24HS-NBFH-S2

Model 36 (45½" L)



CV36LC-5BFH-S1



CV36LS-5BFH-S1



CV36HC-NBFH-S2



CV36HS-NBFH-W2

Model 48 (57½" L)



CV48LC-NBFH-S1



CV48LS-NBFH-W1