


Technical data

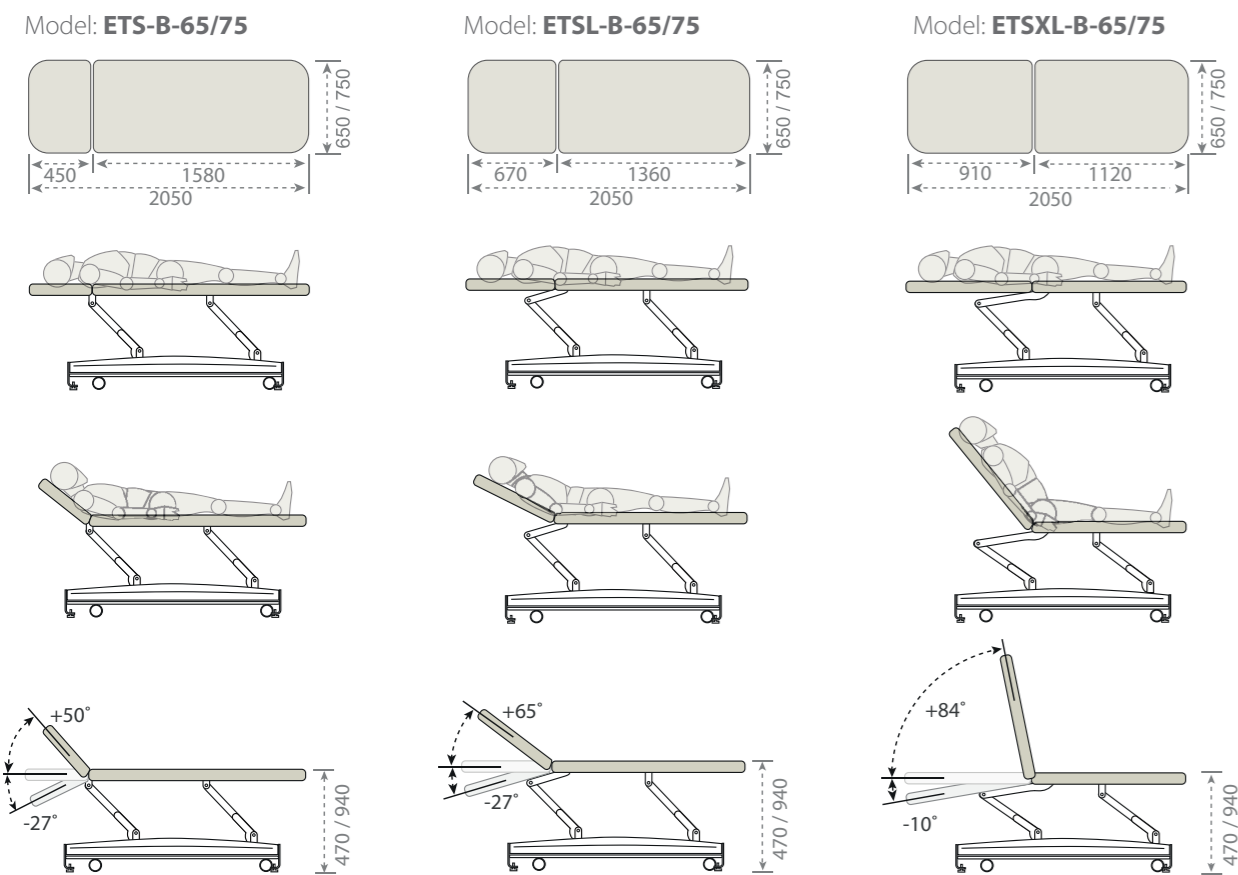


Equipment and accessories

standard ● optional ○

	ETS-B	ETSL-B	ETSXL-B
seamless upholstery (6 cm)	●	●	●
foot switch – SER-FSW	●	●	●
footswitch with safety double click function – SFSW	○	○	○
paper roll holder under head section – NZP1-S, NZP1-SL	○	○	○
fixation rails (pair) – LE-S	○	○	○
infusion stand mounted on fixation rail (pc) – SIF-S	○	○	○
armrest mounted on fixation rail (pc) – NZR	○	○	○
elevated armrest mounted on fixation rail (pc) – NZR1	○	○	○
laminated side table (51x36 cm) mounted on fixation rail – MOL	○	○	○
PVC cover for leg section – PVC6-SB	○	○	○
retractable castors (Ø 75 mm) – SER-RC-SB	●	●	●
antistatic retractable castors (Ø 75 mm, 4 pcs) – ARC-SB	○	○	○
safety switch – EVV-S	○	○	○
plug connector for potential equalization – PPEQ	○	○	○
safe working load 200 kg	●	●	●

Color chart

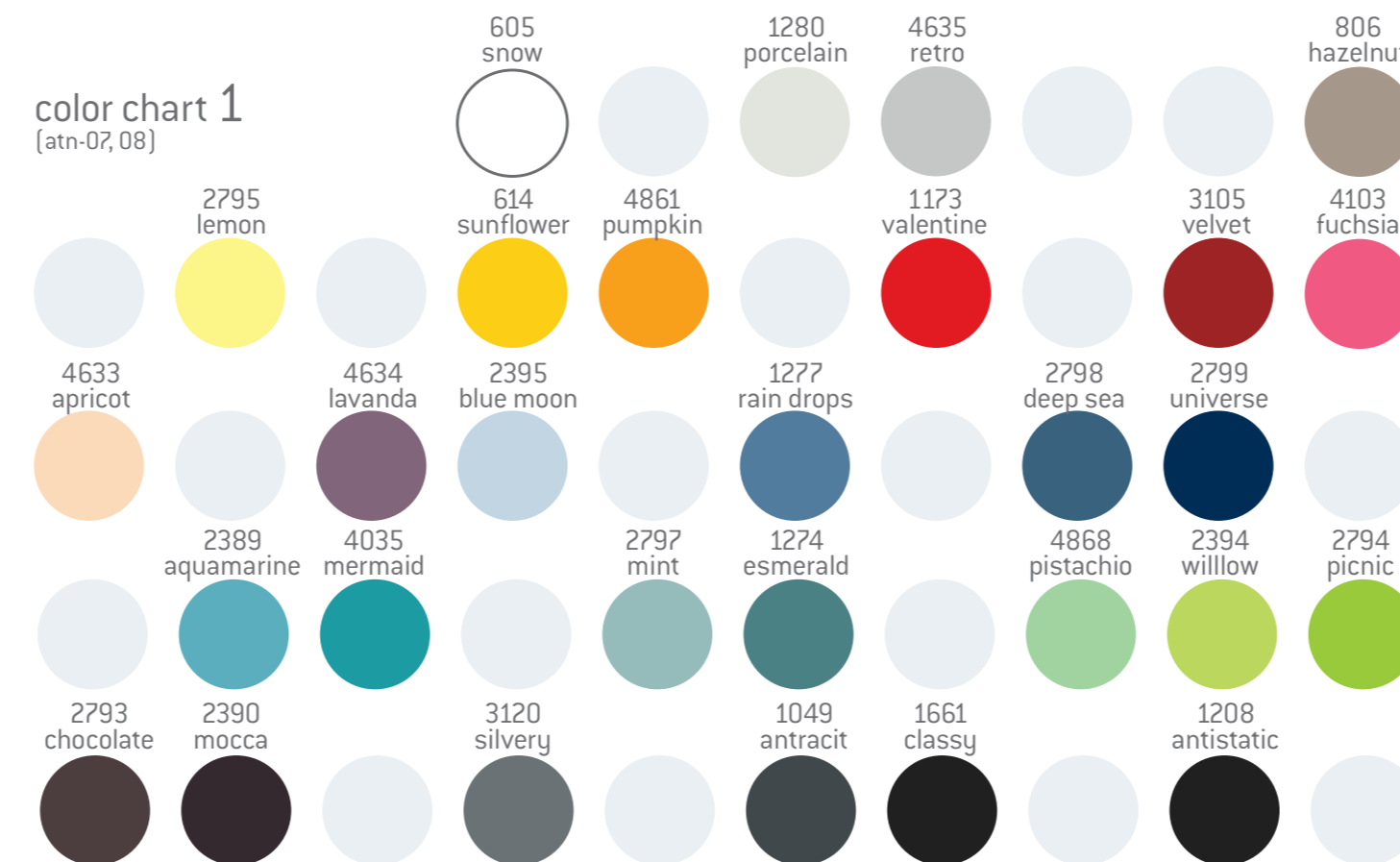

medical examinations



EXAMINATION TABLES S BASIC



Simple adjustment lever design

- easy to press and clean
- stronger due to mechanism integrated in the frame

Upholstery

- seamless 6 cm thick padding for maximum patient comfort
- high quality imitation leather

Retractable castors

- four small castors to easily move the table across the room
- activation lever on both sides of the table

Foot switch

- easy to use
- height adjustment range between 47 and 94 cm

Table base

- protected by durable ABS plastic cover
- attachment is screwless for easier cleaning

Weight capacity

- perfect stability due to standard 200 kg safe working load

Perfect for medical examination

S-line examination tables comprise **THREE BASIC MODELS** distinguished by **flawless quality** and **functionality**. The two-section tables allow **a lying position** and, depending on the head section length, also **a half-sitting** and **sitting position**. Perfectly **STABLE, ELECTRICALLY ADJUSTABLE**, comfortable and with a **wide range of ACCESSORIES**, the tables will fulfill the needs you encounter every day when carrying out tasks **in general or specialized medical practices**.

Safe patient positioning

The **sturdy construction** together with a **200 kg SAFE WORKING LOAD** ensures that examination tables are **perfectly stable** and that the patients can be positioned safely. The tables stand firmly on smartly designed **RUBBER FEET**, but you can easily move them at short distances as they are equipped with **RETRACTABLE CASTORS**. The activation lever is on both sides of the table.

Height adjustment is **electrically operated** by using a **FOOT SWITCH**. It features a simple design and enables you to lift or lower the table to the appropriate height quickly and easily. The height adjustment range is **between 47 and 94 cm**. It allows you to ensure **easy and safe access** for your patients as well as **a comfortable working height** to suit your needs. **Head sections** are **manually adjustable** and supported by gas springs.

Maintenance of high hygiene standards

The examination tables' **simple and clean design** makes it easy to maintain the high hygiene standards required in medical practices. The 6 cm thick **PADDING** is completely **seamless** and the first-class quality **IMITATION LEATHER** is **water- and fire resistant**.

The base is protected by a **COVER** made of durable **ABS PLASTIC**. Its surface is additionally coated with a **scratch resistant coating** and is completely smooth due to **screwless attachment**, which makes the cleaning of the table easier.

Due to the **minimalist design**, the head-section **ADJUSTMENT LEVER** is **easier to press** and **clean** and the mechanism integrated in the frame makes it **stronger**.

medi⁺ HEALTH

Medihealth Oy, Haarasuonkuja 5 90240 Oulu, y-tunnus 3131118-6, www.medihealth.fi

Tilaukset: tilaukset@medihealth.fi Tiedustelut: info@medihealth.fi