

## Description

---

**Stretcher SL** ensures **safety** and **comfort** to the patients and **easy adjustments** and **manageability** to the medical personnel.

It features an **8 cm thick mattress** ensuring perfect patient comfort. It is upholstered in **biocompatible antistatic PU** material resistant to blood, urine, fat and disinfectants. The **welded** and **taped seams** contribute greatly to ensuring **high hygiene standards** as they prevent penetration of liquids.

Hydraulic height adjustment allows a wide **adjustment range between 57 cm and 104 cm**. It provides a very **low access position** for patients and **comfortable height** for the personnel during transport or treatment.

The **sturdy base** protected by the completely smooth **ABS plastic cover** and standard **250 kg safe working load** ensure the patient transport trolley's perfect stability.

Patient safety is guaranteed by a pair of **firm and sturdy safety rails (patent EP 3496694 and US 11110019)** and a possibility of adjustment in **Trendelenburg position** in case of emergency.

Due to **big antistatic castors with central locking system** enabling castors to move freely, in one direction or to be locked, the patient transport trolley is easy to maneuver along the hospital corridors. It features dual height **push handles** at the head section and an additional pair of fixed handles at the leg section. To prevent damage in case of accidental bumping, rubber **bumpers** are at all four corners.



[in](#) [f](#) Bring treatment to a higher level!

**Novak M d.o.o.**

Pod smrekami 5  
SI - 1218 Komenda

E: [info@novak-m.si](mailto:info@novak-m.si)  
T: +386 1 5620117  
F: +386 1 5620118


[www.novak-m.com](http://www.novak-m.com)




---




## Stretcher SL



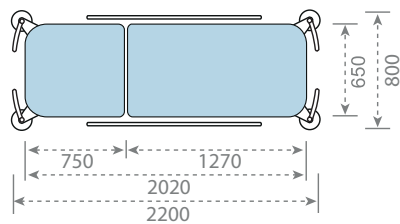


# STRETCHER SL

... a patient transport trolley distinguished by simple and functional design

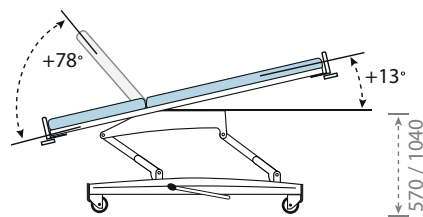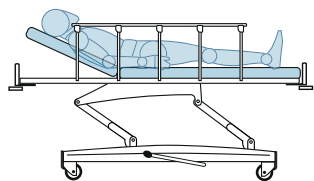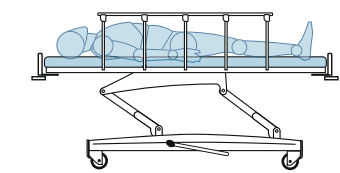

## POSITIONS

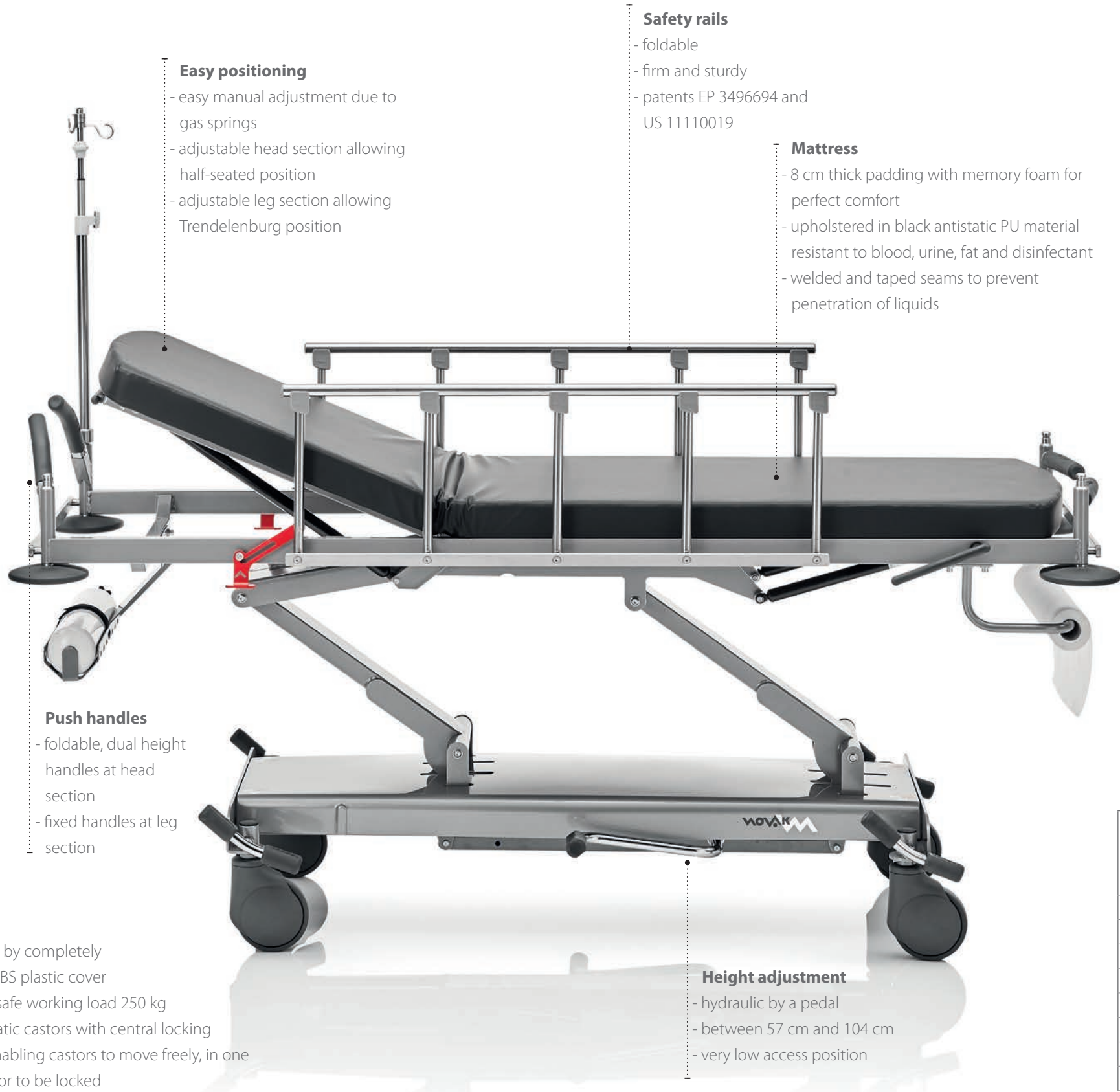

### Easy positioning

- easy manual adjustment due to gas springs
- adjustable head section allowing half-seated position
- adjustable leg section allowing Trendelenburg position

### Safety rails

- foldable
- firm and sturdy
- patents EP 3496694 and US 11110019

### Mattress

- 8 cm thick padding with memory foam for perfect comfort
- upholstered in black antistatic PU material resistant to blood, urine, fat and disinfectant
- welded and taped seams to prevent penetration of liquids

### Bumpers

- on all four corners

### Push handles

- foldable, dual height handles at head section
- fixed handles at leg section

### Base

- protected by completely smooth ABS plastic cover
- standard safe working load 250 kg
- big antistatic castors with central locking system enabling castors to move freely, in one direction or to be locked

### Height adjustment

- hydraulic by a pedal
- between 57 cm and 104 cm
- very low access position



CONFIGURE

Equipment and accessories	STRSL-65
durable antistatic sealed upholstery (8 cm)	●
paper roll holder – NZP2-SL	○
fixation rails at lying, head or leg section – LE2-S, LE-X-S, LE-X-Z	○
infusion stand – SIF-X	○
oxygen bottle holder – O2-S (45 cm)	○
oxygen bottle holder – O2-S-L (55 cm)	○
antistatic castors with central locking system (Ø 150 mm)	●
safe working load 250 kg	●

standard ● optional ○

# medi<sup>+</sup> HEALTH

Medihealth Oy, Haarasuonkuja 5 90240 Oulu, y-tunnus 3131118-6, [www.medihealth.fi](http://www.medihealth.fi)

Tilaukset: [tilaukset@medihealth.fi](mailto:tilaukset@medihealth.fi) Tiedustelut: [info@medihealth.fi](mailto:info@medihealth.fi)