


Technical data

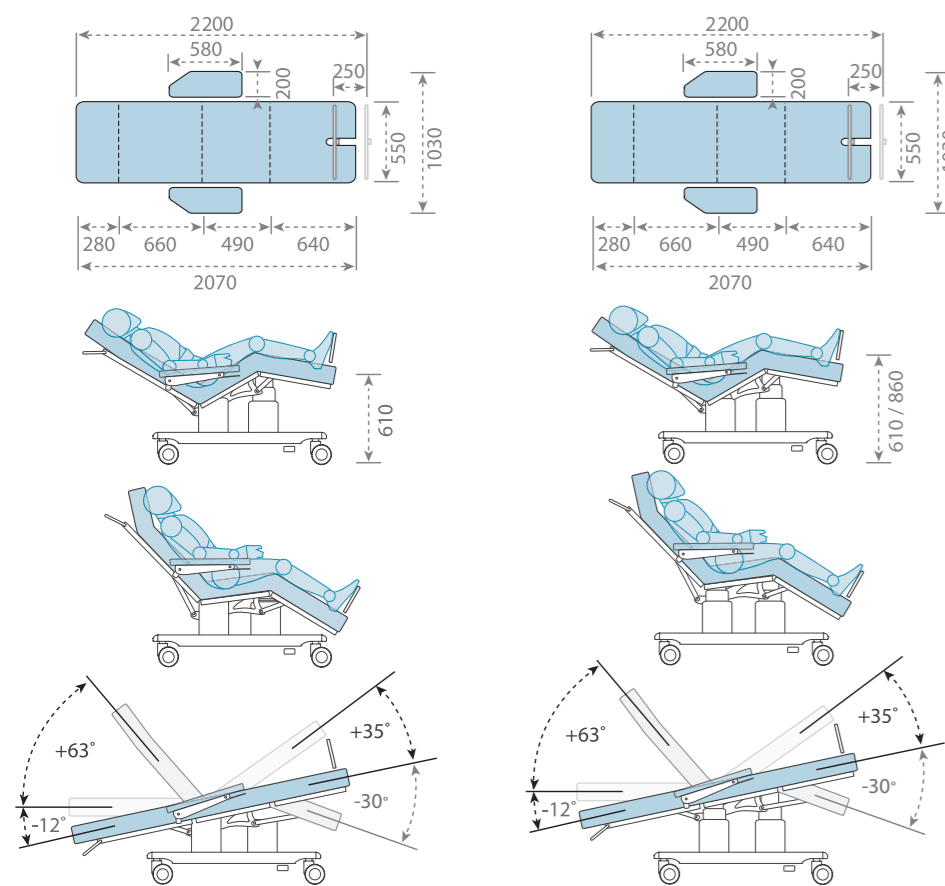


Equipment and accessories

	TCZ-EX	TCZ-LUX
height of the chair - upholstery included	61 cm	61-86 cm
height of the chair - upholstery excluded	47 cm	47-72 cm
safe working load	250 kg	250 kg
push handles	●	●
paper roll holder – NZP1-Z	○	○
adjustable footrest – FR-Z	○	○
cushion for footrest – BZNI-Z	○	○
infusion stand holder, left – SIFHI-Z-L	○	○
infusion stand holder, right – SIFHI-Z-R	○	○
infusion stand – SIF1-Z	○	○
fixation rail at sitting section (pair) – LE1-Z	○	○
armrest side protection (pair) – HRP-Z	○	○
PVC cover for leg section – PVC1-Z	○	○
PVC cover for leg section with footrest – PVC1-FR-Z	○	○
PVC cover for armrest – PVC3-Z	○	○
castors with central locking system (Ø 100 mm) – SER-CK1-Z	●	●
antistatic castors with central locking system (Ø 100 mm) – AC1-CK-Z	○	○
large antistatic castors with central locking system (Ø 125 mm) – AC2-CK-Z	○	○
plug connector for potential equalization – PPEQ	○	○
battery backup – BB-Z	○	○

standard ● optional ○

Color chart PU

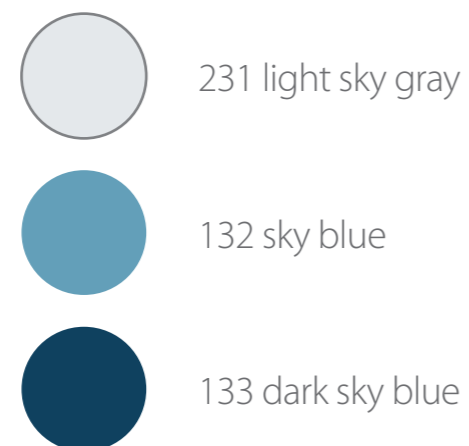

Therapy chair ZeroG

Where time flies!



push handles



manually adjustable head section

extremely soft padding for optimal support (14 cm)

simply adjustable armrests

castors with central locking system (Ø 100 mm)

Therapy chair **ZeroG** exclusive
Model: TCZ-EX

»The hand switch is easy to operate and the illustrations are clear and understandable. When I put the chair in ZeroG position, I didn't move until the end of the therapy. Due to the exceptionally soft and comfortable upholstery, I was calm and without back pains. It felt like I was sitting on a soft cloud.«
(Maja, patient at dialysis department)



easy adjustable footrest (optional)



infusion stand (optional)

paper roll holder (optional)

infusion stand holder, left (optional)

fixation rail (optional)

electric height adjustment

PVC cover for leg section with footrest (optional)

antistatic castors with central locking system (Ø 100 mm) (optional)

Therapy chair **ZeroG** luxury
Model: TCZ-LUX

»Therapy chair ZeroG is a luxury for the patients, as they can set it into desirable position all on their own. Height adjustment is luxurious especially for the medical staff, because we can take better care of our back.«
(Mojca, nurse at rheumatology department)

Comfort

The **soft padding** and setting the therapy chair into **Zero Gravity position** guarantee first-class comfort for patients during long-lasting therapies. The upholstery is made of **high-quality foam** that provides **full body support** along the lying surface and prevents **overheating** and **sweating**. By pressing the **“cloud” button** on the hand switch, the patient can set the chair into **ZeroG** position. This offers additional **pressure relief for the body and joints** during longer therapies.

Easy adjustment

A wide range of smart and simple features save handling time for medical staff, as patients can **adjust the chair independently** into their favourable position. By pressing a single button on a hand switch, the patient can set the chair into **ZeroG position, sitting position** or **completely flat position**. In case of emergency you can quickly set the chair into **Trendelenburg** position. All chair sections can be adjusted **simultaneously**. The chair can be additionally equipped with **footrest** that supports patients in case of cramps. By pressing and releasing the built-in button they can quickly adjust the footrest themselves. For even greater comfort, the patient can easily adjust the **armrests** by pressing the lever. Besides offering support during therapies, strong and stable armrests also help patients with sitting and getting up from the chair.

Simple cleaning and maintenance

Therapy chair ZeroG has a **small footprint** so it takes very little space in a clinic. Simply designed lying surface is smooth and without folds for easier and faster cleaning. It is made of thin, **extra strong PU material**, which is resistant to blood, urine, oil, fat and disinfectants to **ensure proper hygiene** for your patients.

medi⁺ HEALTH

Medihealth Oy, Haarasuonkuja 5 90240 Oulu, y-tunnus 3131118-6, www.medihealth.fi

Tilaukset: tilaukset@medihealth.fi Tiedustelut: info@medihealth.fi